

Stay protected
24/7



The **SkyHelp** is a mobile medical alert system. It is designed to keep you safe, at home or on the go. There is no base station or landline needed to use your **SkyHelp**. The **SkyHelp** has a microphone and speaker to enable easy two-way communication with one of our highly-trained emergency operators. The **SkyHelp** has 3G wireless coverage so you can take it anywhere with cellular coverage and get help when you need it. With SkyHelp's GPS Tracking System we can find you when you cannot speak.

Lightweight - 1.5 ounces - comfortable to carry or wear.

Durable - Shatter-resistant with rugged construction.

Small In Size - 2.5" x 1.8" x .67 - compact yet powerful.

Loud - Crystal clear audio quality - No need to hold to the ear or mouth to listen or speak.

One-Touch Operation and Braille Button - Press the button 2 seconds and the call is placed. Anti-cancel prevents cancelling the call if the button is pressed again. 'Help' in braille is written on the button as well for extra protection.

Waterproof - Wear it in the bath or shower where most accidents occur. *IPX7 rating states product is protected against high pressure water jets from any direction. Limited ingress permitted (ex. shower). Protected against short periods of immersion in water (3.3 feet of water for up to 30 minutes - ex. bath).

Wireless Network - Take it wherever you go on a nationwide 3G wireless network

Voice Prompts - Lets you know the status of your call. Ex. Your emergency call is now being dialed, please stand by.



Direct GPS Location - Sends your position directly to the monitoring center and/or a loved one.


Fall Detection - Notifies the monitoring center/loved ones if you fall.

Geo-Fencing - Notifies operator/friends/family if you are near and/or step out of your pre-programmed 'safe zone'.

Extended Battery Life - 72 Hour Use - For maximum use between charges.

Battery Status Indicator/Fully Supervised Battery - Visual and audio indicator of battery status. Battery indicator lights (and cellular strength indicator light) surround the call button and will let you know when the battery or cellular signal strength is low. The pendant will also announce that the battery is low, and the unit will notify the monitoring center as well. The monitoring center will then report the low battery status back to you.

Protected Phone Number - Only the monitoring center can call you on your protected phone number, and this will ring your SkyHelp. This will help to locate your pendant in the event of misplacement.

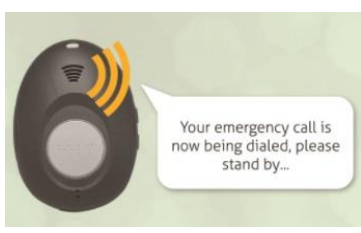

- 1 SkyHelp Pendant
- 1 Docking/Charging Station
- 1 Power Cord
- 1 Power Cord Adapter
- 1 12 Volt Adapter
- 1 Belt Clip
- 1 Lanyard





ORDER FORM

Assistive Technology Services
"Helping You Live Independently"

8023 Franklin Road
 Murfreesboro, TN, 37128
 (615) 562-0043
 Fax (831) 308-9671

Billing Address Name _____ Address _____ City _____ State _____ Phone Number _____	Check if same as Billing [<input type="checkbox"/>]	Ship To Address Name _____ Address _____ City _____ State _____ Phone Number _____
---	--	---

CUSTOMER EMAIL ADDRESS	PAYMENT TERMS	DATE
	Due on receipt	

QTY	ITEM #	DESCRIPTION	UNIT PRICE	SHIPPING	UNIT TOTAL
	atsSkyHelp	SkyHelp Monthly Plan	\$179.99	Included	
		Includes Setup Charge, Shipping, Activation, SkyHelp, Charging Cord, Charging Station, Power Adapter, Lanyard, Belt Clip and First Month Service	\$39.99/mo		
	atsSkyHelp	SkyHelp Annual Plan	\$179.99	Included	
		Includes Setup Charge, Shipping, Activation, SkyHelp, Charging Cord, Charging Station, Power Adapter, Lanyard, Belt Clip and One Year Service	\$449.99/yr		
				Subtotal	
		Only Applicable If Shipping To Tennessee		TN Sales Tax 9.75%	
		*Prices are Subject to Change without Notice		Shipping	
				Total	

If Paying By Check Please Make Checks Payable To: Assistive Technology Services

Credit Card Type	Credit Card Number	Expiration Month	Expiration Year	Security Code From Back Of Card
Visa [<input type="checkbox"/>] MC [<input type="checkbox"/>] AmEx [<input type="checkbox"/>] Disc [<input type="checkbox"/>]				
Name On Card	Email Address			

"May You Know The Glory Of God And May He Bless You"